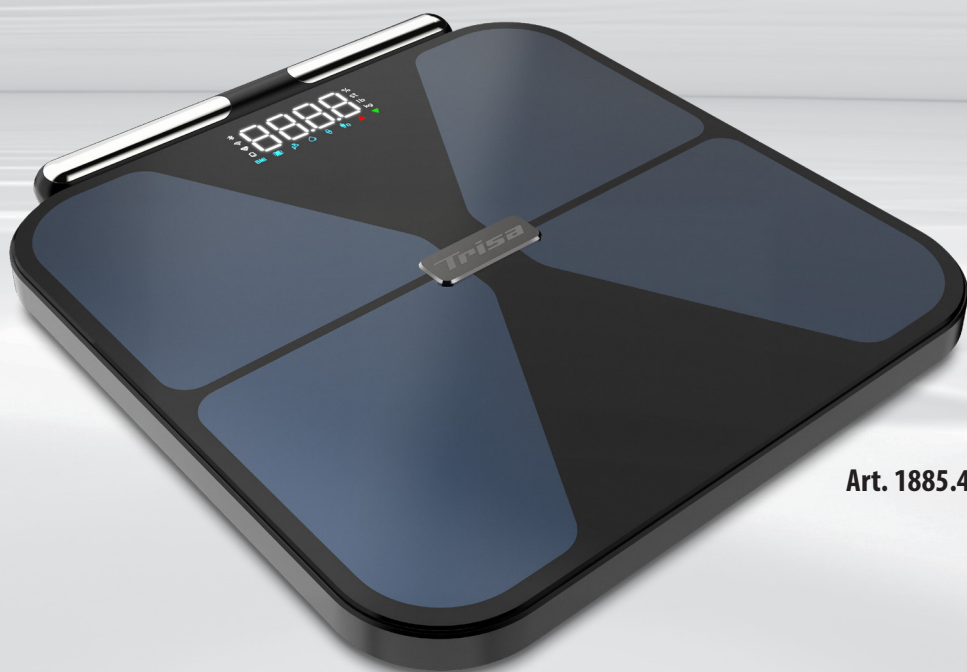
Art. 1885.42

Personenwaage
Bathroom scales
BioMetrix Smart

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EN – Safety instructions

Read the safety instructions carefully before using the device and keep them for future reference.

- This scale is not suitable for people with pacemakers or other medical implants. The indication of body fat may be inaccurate in the event of diabetes or other medical/physical conditions. The same applies to highly trained sportsmen.
- This device is not suitable for pregnant women!
- Any treatment or diet for overweight or underweight requires qualified advice from appropriate specialists (doctor, dietician). The values determined with the scale can be of help.
- Risk of tilting! Never place the scales on an uneven surface. Never step onto just one corner of the scale. While weighing, never stand on one side or on the edge of the scale.
- Danger of slipping! Never step onto the scale with wet feet. Never step onto the scale while wearing socks.
- Use the device only for its intended purpose, as directed by the user instructions. In the event of misuse, all guarantees become null and void.
- The scale is manufactured for home use. It is not intended for commercial use in hospitals or other medical facilities.
- Do not use the device if it does not work properly, if it has been dropped, dropped into water or damaged.
- The scale has a measuring range up to 180 kg/396 lbs. Do not overload the scale.
- Place the scale on a firm level surface. Soft, uneven surfaces are unsuitable for measurement and lead to incorrect results.
- Place the scale in a location where neither extreme temperatures nor extreme humidity are present.
- Keep the scale away from water.
- Please handle the scale with care. Avoid shocks and vibrations of the scale. Do not drop it.
- Stand carefully on the scale. Do not hop or jump on the tread surface, as this may damage the weighing mechanism.
- Do not attempt to disassemble the device, otherwise the warranty will become void. This device does not contain any

parts that can be serviced or replaced by the user.

- In the event of malfunctions, do not proceed with any unauthorised repair of the device. Otherwise, this will render any warranty claims void. Have repairs performed only by an authorised service centre.
- After taking the measurement, carefully lower the handle and do not release it beforehand.
- WARNING: Make sure that children do not play with the packaging films. Choking hazard!

Battery safety instructions

- Do not take the batteries apart!
- If required, clean the battery and device contacts before insertion!
- Removed spent batteries from the device immediately!
- Increased risk of leakage, avoid contact with skin, eyes and mucous membranes!
- In the event of contact with battery acid, immediately rinse the affected area with plenty of clean water and consult a doctor immediately!
- If a battery is swallowed, consult a doctor immediately!
- Insert the batteries correctly, paying attention to the polarity!
- Keep the battery compartment closed properly!
- If the device is not used for an extended period of time, remove the batteries!
- Keep the batteries out of the reach of children!
- Do not re-charge the batteries! There is a risk of explosion!
- Do not short-circuit the device! There is a risk of explosion!
- Do not throw into a fire! There is a risk of explosion!
- Store unused batteries in the packaging and not near metal objects to avoid a short circuit!

Scope of delivery

First, check whether the device is complete and is not damaged. If in doubt, do not operate the device and contact your service centre.

Disposal

Render worn-out appliances unusable. Unplug and cut off cord. Dispose of electrical appliances at collection points or in store. Do not dispose of with household waste. Follow local regulations.

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Safety Precautions

Read these safety instructions carefully before use.

Product Introduction

This product features a LED display and high-precision sensors for weight measurement. It supports automatic step-on activation and can measure up to 56 body composition metrics.

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
Precautions

- Place the scale on a flat, hard, non-slippery surface.
- Do not use if pregnant or with implanted medical devices (pacemaker, artificial lung, wearable ECG).
- Note: To measure body composition beyond weight, step barefoot onto the scale.
- Due to the natural growth process, measurement results for persons under the age of 18 should be considered reference values only.

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Download **Unique Health** App via app store or QR code.

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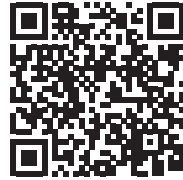Follow the app instructions to register
and log in.

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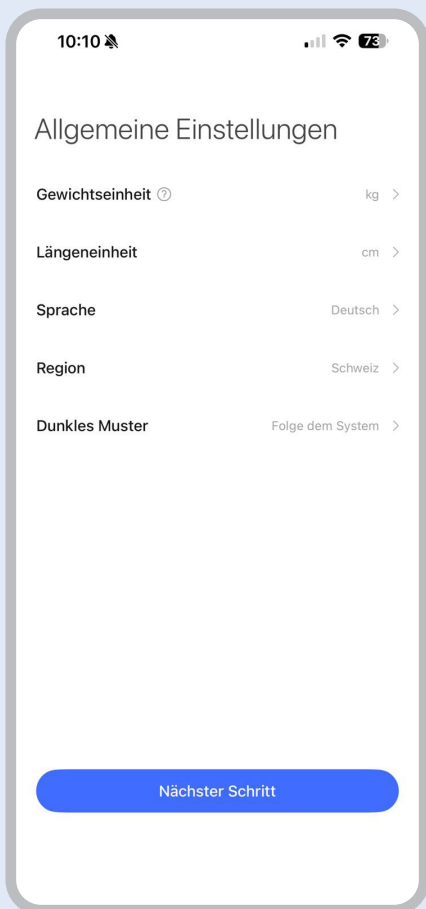
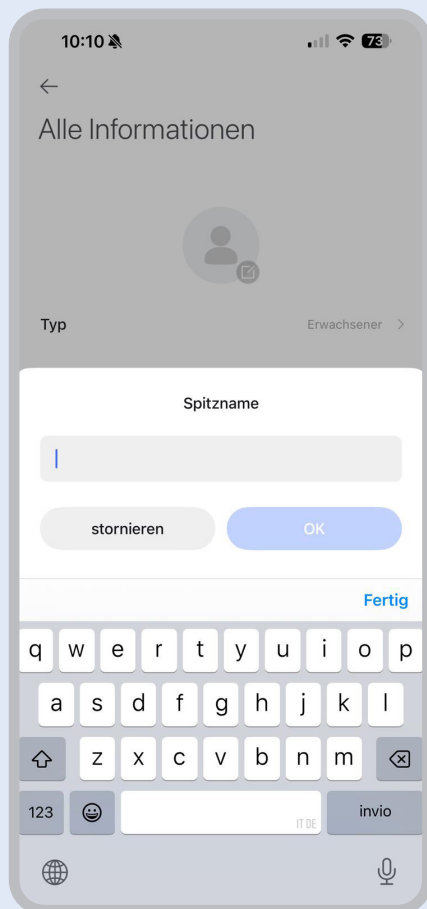
Request the verification code by tapping “Get verification code”.

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When the 60-second timer appears, the code will be sent by email.

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
Check your email, copy the code, and paste it into the app.

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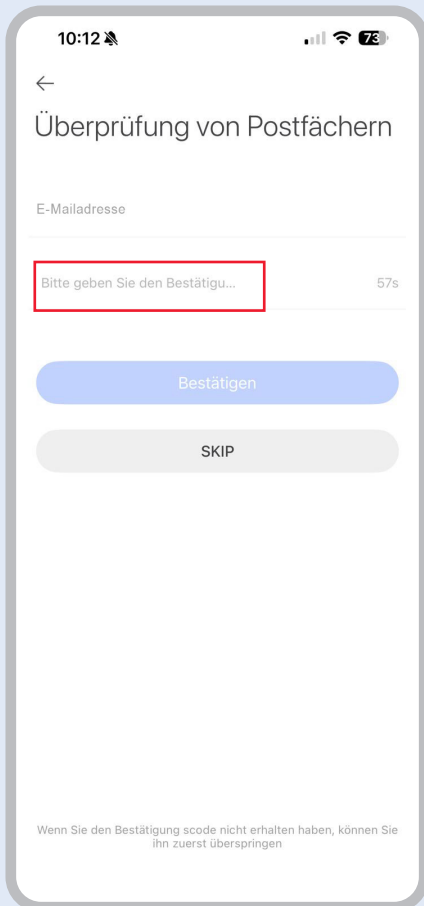
Login completed successfully.
The app is now ready to use.

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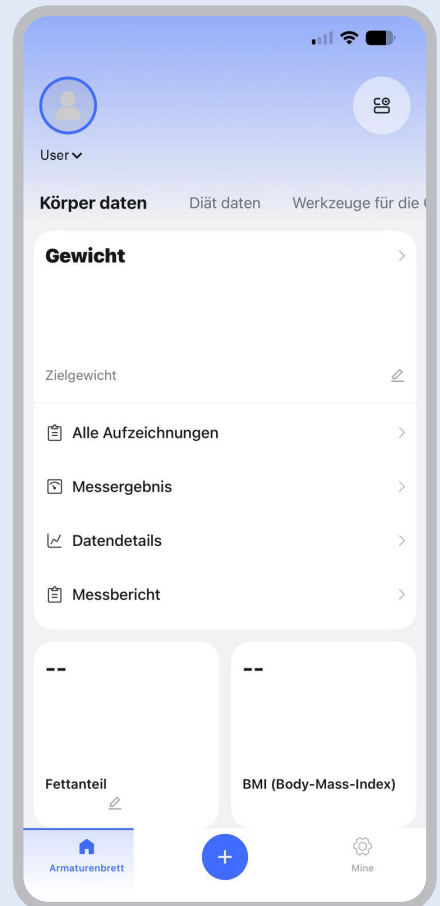
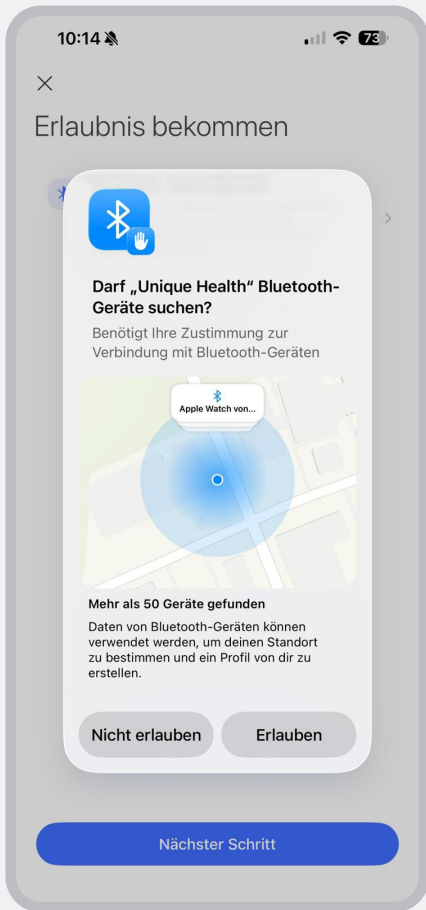
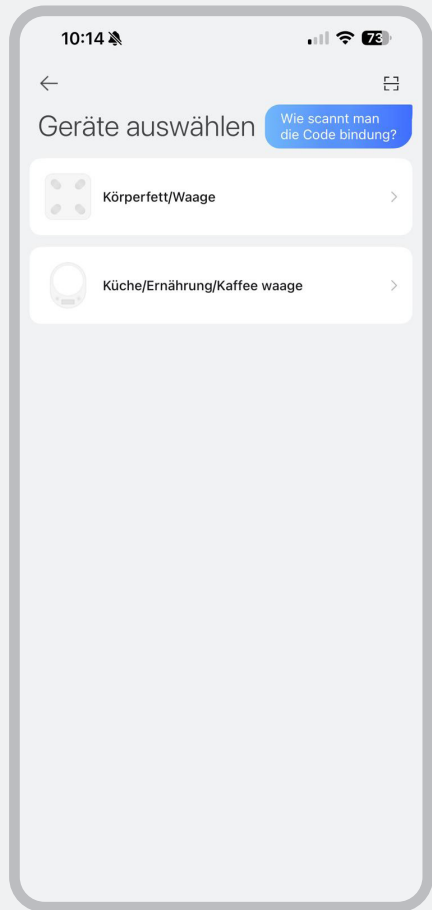
Step onto the scale briefly to activate it.

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
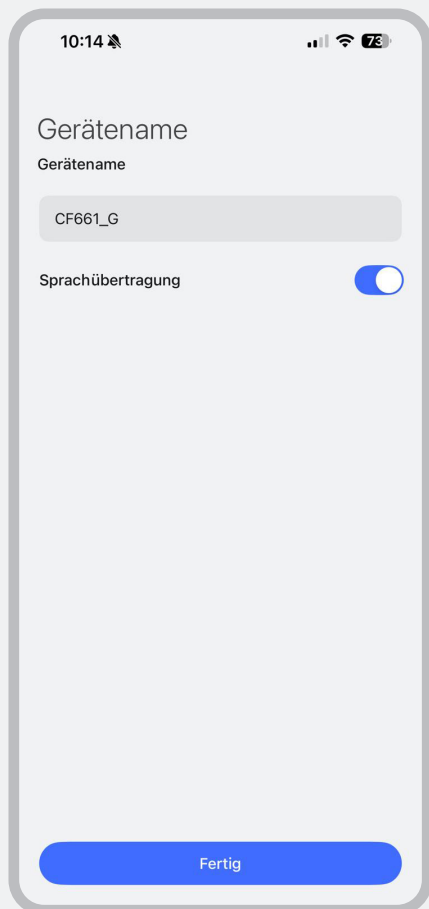
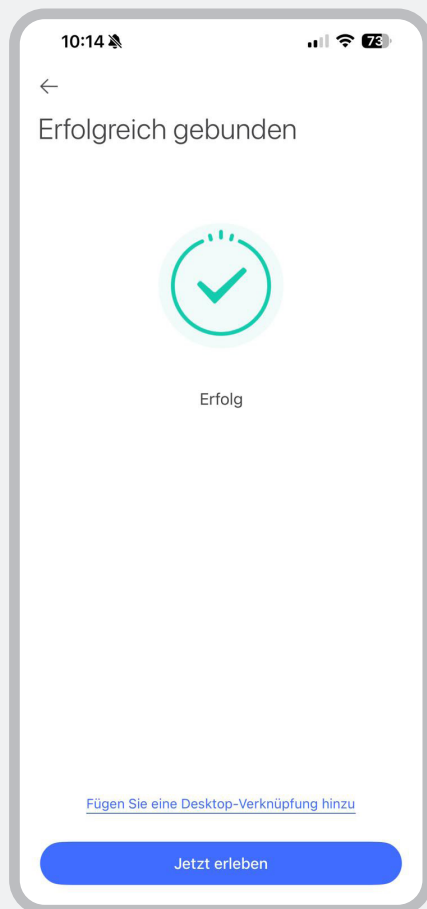
Connect to the network and continue with the configuration.

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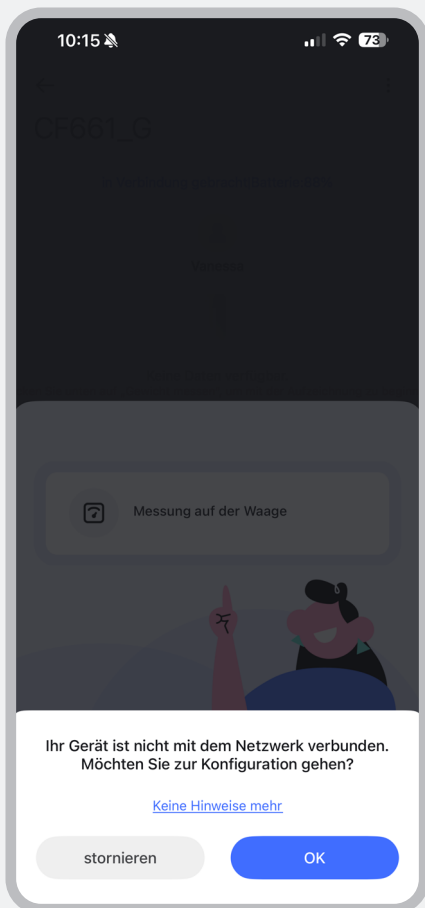
Configure the network and connect to your own Wi-Fi network.

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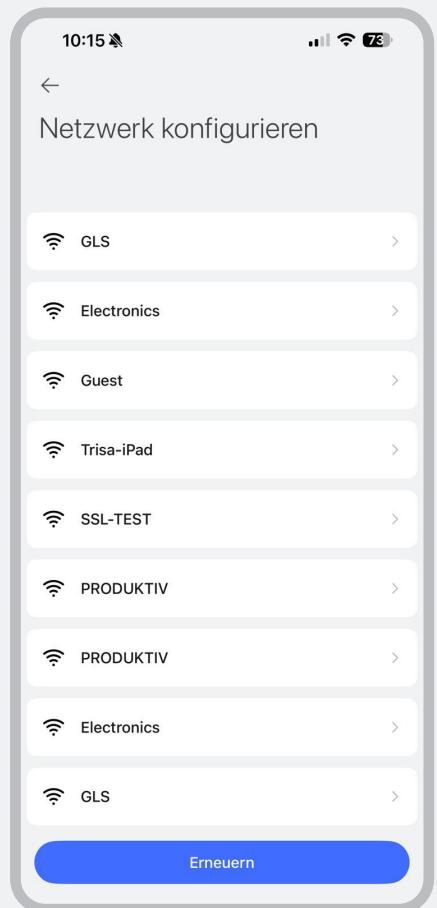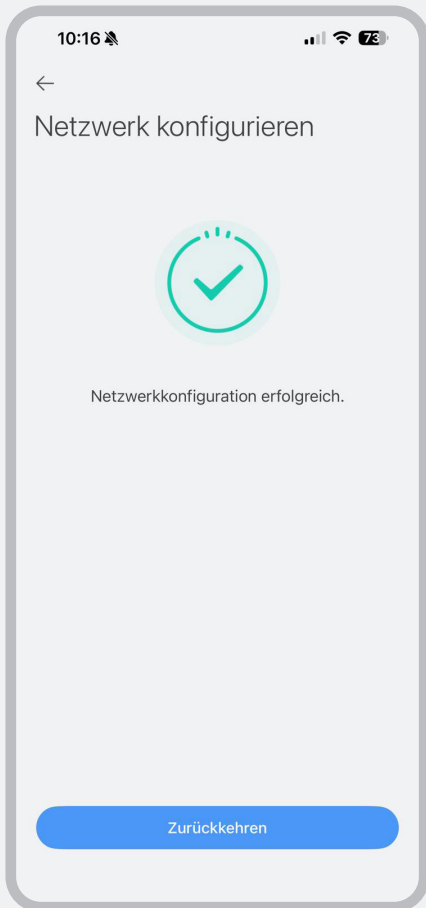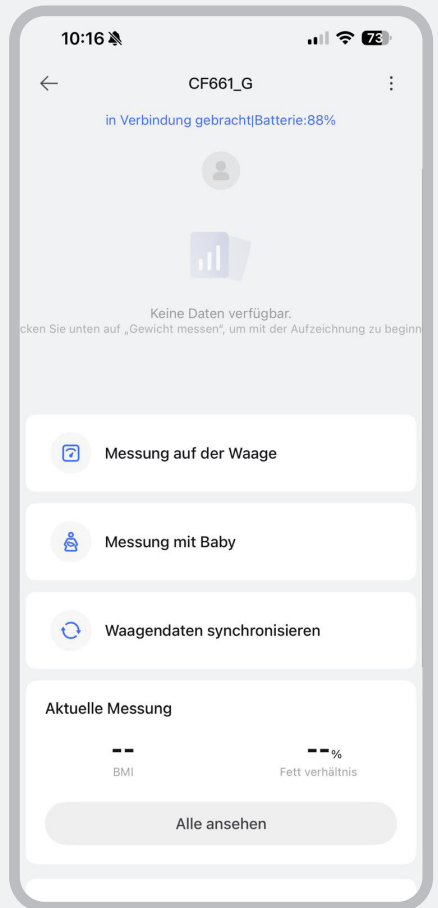
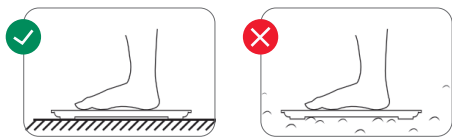
Place the scale on a hard, flat surface. Do not place the scale on carpets or soft flooring.

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Place both feet evenly and centrally on the scale.

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
Pull the handle upward, hold it firmly, and extend the arms, as shown in the illustration.


After the measurement, lower the handle slowly and do not release it suddenly.

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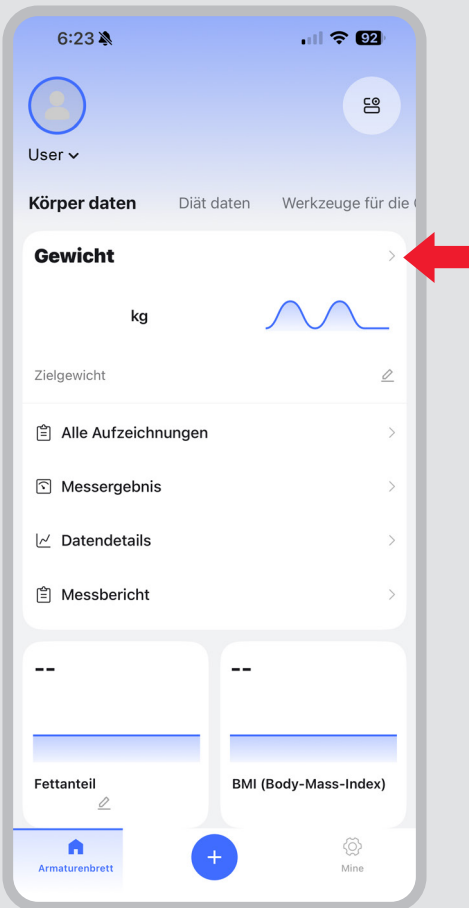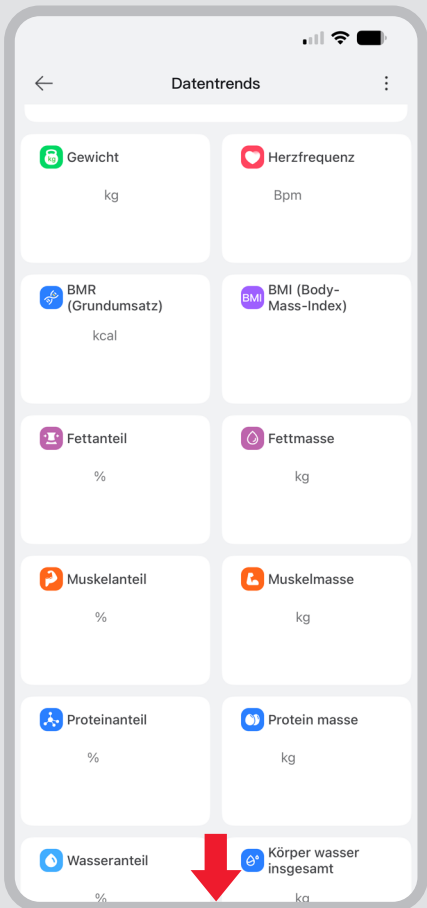
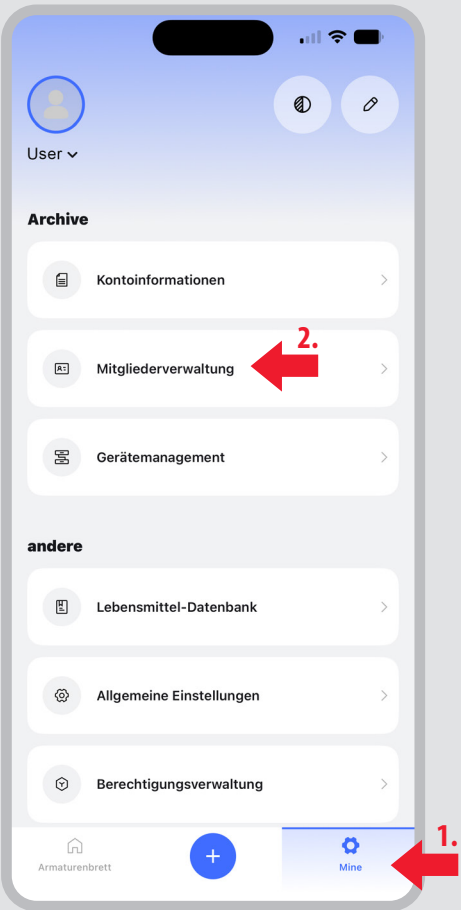
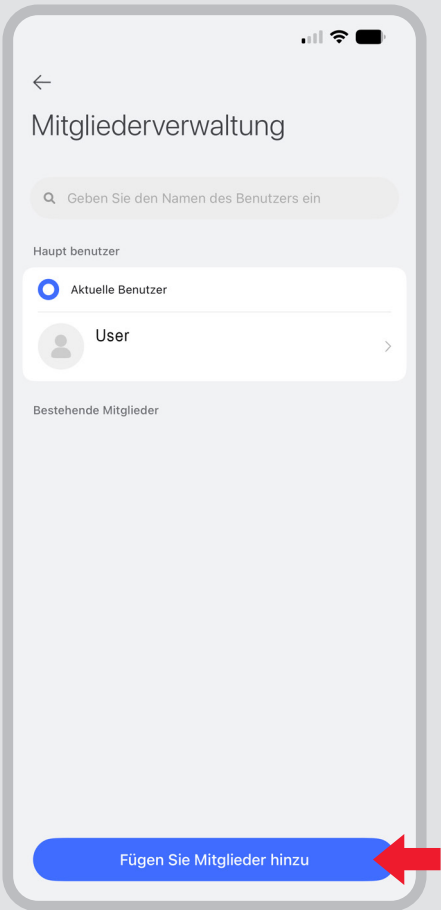
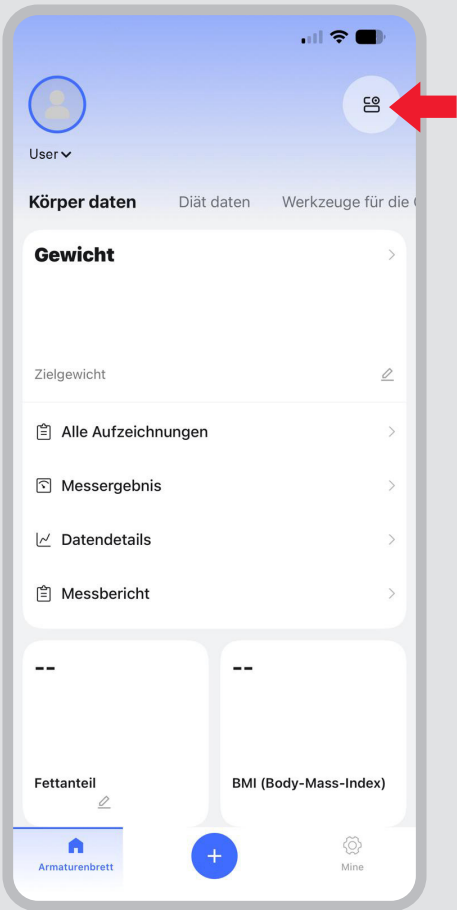
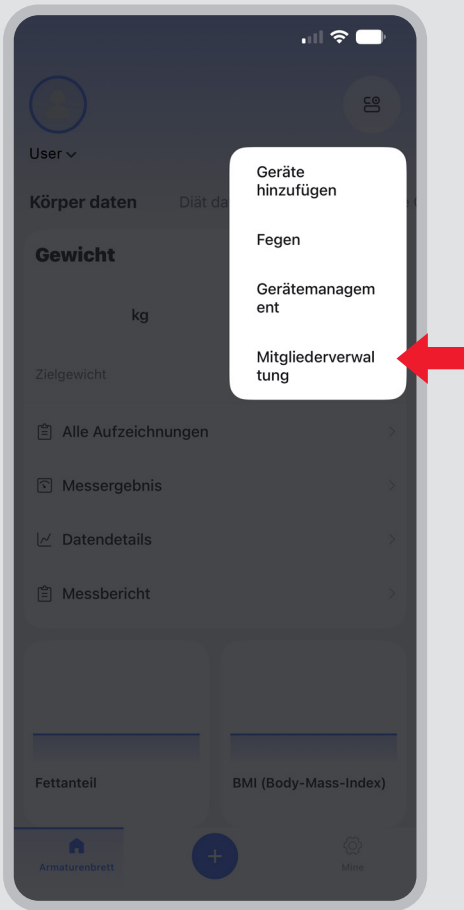
Insert batteries: Open the battery compartment on the bottom and insert the batteries.

Remove batteries if not used for a long time. Use suitable batteries only and dispose of them properly.

Anzeige | Display

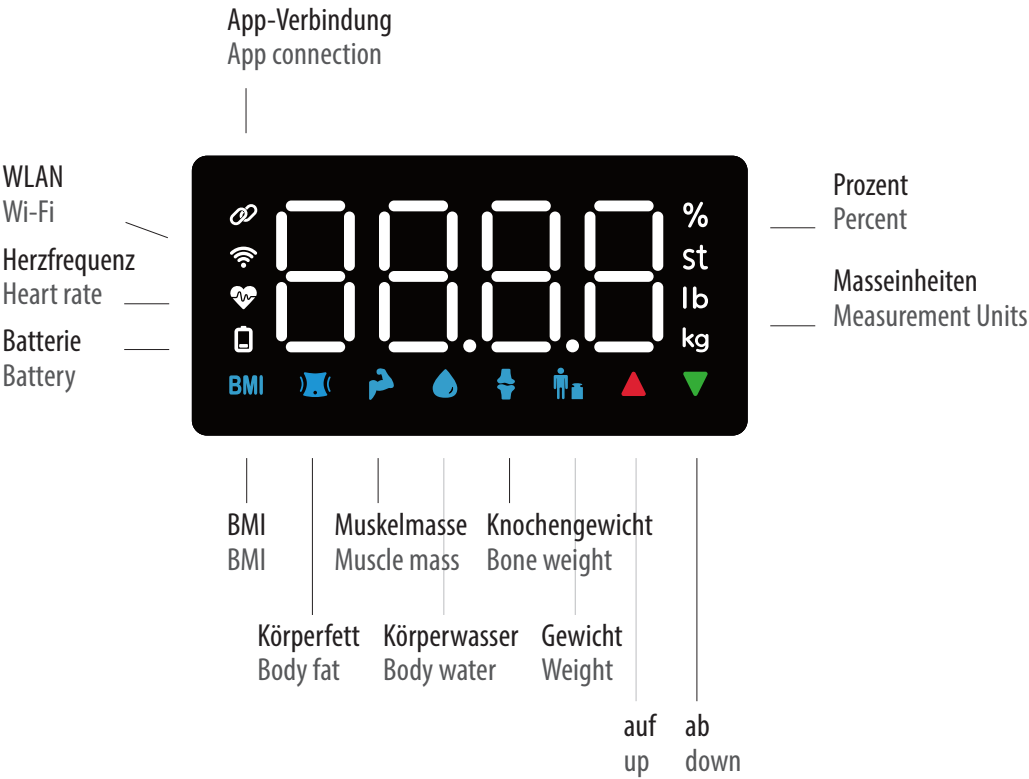CH

Trisa Electronics AG
Kantonsstrasse 121
CH-6234 Triengen
info@trisaelectronics.ch
+41 41 933 00 30

AT

SCD Handels GmbH
Unterhaus 33
A-2851 Krumbach
service@trisaelectronics.at
+43 (2647) 430470

DE

Ernst Lorch KG
Bildstockstrasse 9
D-72458 Albstadt
HG@lorch-gruppe.com
+49 (7431) 124 200

BG

Trisa GmbH
43A, Gorski Patnik Strasse
BG-1421 Sofia

Service & Support
www.trisaelectronics.ch